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Wundversorgung

EINFACH.

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Referenzen Safetac® Technologie:

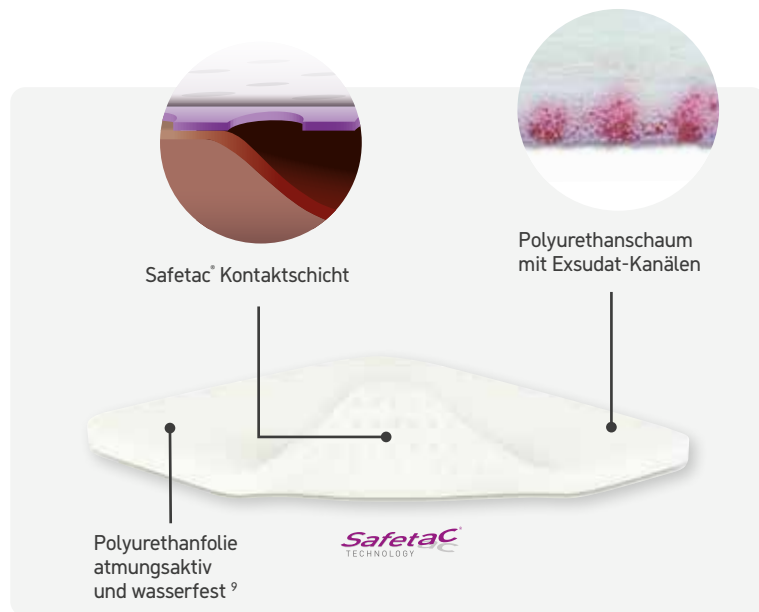
1. Van Overschelde, P. et al. A randomised controlled trial comparing two wound dressings used after elective hip and knee arthroplasty. Poster presentation at 5th Congress of the WUWHs, Florence, Italy, 2016. 2. Silverstein P. et al. An open, parallel, randomized, comparative, multicenter study to evaluate the cost-effectiveness, performance, tolerance, and safety of a silver-containing soft silicone foam. Journal of Burn Care and Research, 2011. 3. Gee Kee E.L. et al. Randomized controlled trial of three burns dressings for partial thickness burns in children. Burns, 2014. 4. David F. et al. A randomised, controlled, non-inferiority trial comparing the performance of a soft silicone-coated wound contact layer (Mepitel One) with a lipidocolloid wound contact layer (UrgoTul) in the treatment of acute wounds. International Wound Journal, 2017. 5. Patton M.L. et al. An open, prospective, randomized pilot investigation evaluating pain with the use of a soft silicone wound contact layer vs. a traditional dressing and staples on split thickness skin grafts as a primary dressing. Journal of burn care & research, 2013. 6. Bredow J. et al. Evaluation of Absorbent Versus Conventional Wound Dressing. A Randomized Controlled Study in Orthopedic Surgery. Deutsche Ärzteblatt International, 2018. 7. Meaume S. et al. A study to compare a new self-adherent soft silicone dressing with a self-adherent polymer dressing in stage II pressure ulcers. Ostomy Wound Management, 2003. 8. Gotschall C.S. et al. Prospective, randomized study of the efficacy of Mepitel on children with partial-thickness scalds. Journal of Burn Care & Rehabilitation, 1998.

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Tubifast®
WITH
2-WAY STRETCH®
TECHNOLOGY

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Tubifast®
WITH
2-WAY STRETCH
TECHNOLOGY

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Referenzen:

1. Franks PJ, Moody M, Moffatt CJ, Hiskett G, Gatto P, Davies C, et al. Randomized trial of two foam dressings in the management of chronic venous ulceration. Wound repair and regeneration: official publication of the Wound Healing Society [and] the European Tissue Repair Society. 2007;15(2):197-202. 2. Mölnlycke Health Care. Data on file. 2016. 3. White R. A multinational survey of the assessment of pain when removing dressings. Wounds UK 2008; 4(1):14-22. 4. Schumann, H., Beljan, G., Hoping, D., Bruckner-Tuderman, L. Atraumatic dressings in fragile skin conditions: use of the soft silicone dressing (Mepilex) in hereditary and acquired bullous skin disease. Poster presentation. EWMA. 2005. 5. Upton, D., Solowiej, K. The impact of atraumatic vs conventional dressings on pain and stress. Journal of Wound Care 2012;21(5):209-215. 6. Mölnlycke Health Care. Data on file. 2017. 7. Weaver, G., Crawford, K. The use of Mepilex on dehiscent amputation wounds. Wounds UK 2007;3(1):70-74.

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Referenzen:

1. Meuleneire, F., and Fostier, A. Local treatment of heel pressure ulcers with a silicone foam dressing. Poster presentation, WUWHS, 2008. 2. Perez YP, Carmona JA, Perez IL, Garcia CM. Prevention and treatment of radiodermatitis using a non-adhesive foam dressing. Journal of Wound Care 2011;20 (3):130-135. 3. Diggelmann KV, Zytkowicz AE, Tuaine JM, Bennett NC, Kelly LE, Herst PM. Mepilex Lite dressings for the management of radiation-induced erythema: a systematic inpatient controlled clinical trial. British Journal of Radiology, 2010;83(995):971-8. 4. Paterson DB, Poonam, P., Bennett, NC, Peszynski, RI, Van Beekhuizen, MJ, Jasperse, M, Herst, PM. Randomized intra-patient controlled trial of Mepilex Lite dressings versus aqueous cream in managing radiation-induced skin reactions post-mastectomy. Journal of Cancer Science and Therapy, 2012;4 (11):347-56. 5. 24. Boesch RP, Myers C, Garrett T et al. Pediatrics, 2012; 129: e792-e797. 2. 2018. 6. Zheng XP, Huang GY, Chang F, Qian MY, Xia ZF, Xiao SC. Curative effect of soft silicone dressing combined with calcium alginate dressing in treating skin graft donor sites of burned patients. Academic Journal of Second Military Medical University, 2016;37(11):1321-4. 7. Schumann, H., Beljan, G., Hoping, D., Bruckner-Tuderman, L. Atraumatic dressings in fragile skin conditions: use of the soft silicone dressing (Mepilex) in hereditary and acquired bullous skin disease. Poster presentation: EWMA, 2005. 8. White R. A multinational survey of the assessment of pain when removing dressings. Wounds UK 2008; 4(1):14-22. 9. Upton, D., Solowiej, K. The impact of atraumatic vs conventional dressings on pain and stress. Journal of Wound Care 2012 21(5):209-215. 10. Eytier C, Gazeau E, Beneteau G, Verfaillie G. Convenience and tolerance of the combination of a soft silicone foam dressing and a two-way stretch tubular bandage in the management of local wounds. Journal des plaies et cicatrisations 2013;18(88):38-44. 11. Mölnlycke Health Care. Data on file. 2006. 12. Zhang Y, Xing SZ. Treatment of Diabetic Foot Ulcers using Mepilex Lite Dressings: A Pilot Study. Experimental and Clinical Endocrinology and Diabetes: official journal, German Society of Endocrinology [and] German Diabetes Association. 2014;122(4):227-30. 13. Mölnlycke Health Care. Data on file. 2016. 14. Khramilin V. Mepilex Lite/EM in the treatment of diabetic foot ulcer. Poster presentation, EWMA, 2006. 15. Mölnlycke Health Care. Data on file. 2016. 16. Mölnlycke Health Care. Data on file. 2018.

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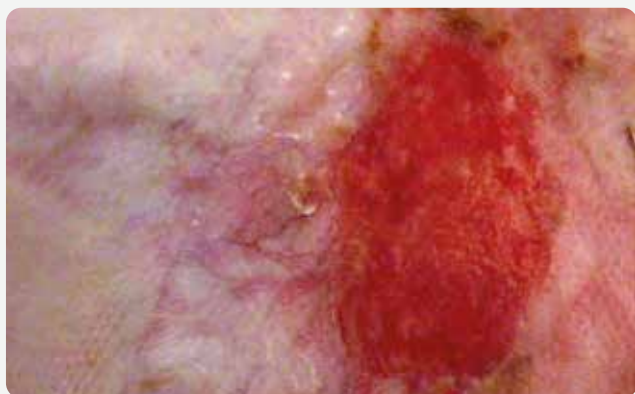
Referenzen:

1. Product Manual – Conformability PD-528870 Data on file 2. External data – FEM simulations Mx-B Flex vs Allevyn Life PD-529747 3. Product Manual – Fluid Handling capacity PD-527642. Data on file 4. Product Manual – Viscous fluid PD-528871. Data on file 5. Product Manual – Mepilex Border Flex Bacteria encapsulation PD-537072. Data on file 6. Product – Viral penetration. Report no. PD-535090. 7 Karlsson, C., Bianchet, A., Blom, K., Almegren, E., Wellner, E., Hamberg, K. Bacterial trapping of a newly developed all-in-one soft silicone foam dressing. Poster presentation at European Wound Management Association conference, Krakow, Poland, 2018 7. Santamaria, N., Gerditz, M., Kapp, S., Wilson, L., Gefen, A. A randomised controlled trial of the clinical effectiveness of multi-layer silicone foam dressings for the prevention of pressure injuries in high-risk aged care residents: The Border III Trial. International Wound Journal 2018 DOI: 10.1111/iwj.12891. 8. Kalowes, P., Messina, V., Li, M. Five-layered soft silicone foam dressing to prevent pressure ulcers in the intensive care unit. American Journal of Critical Care 2016;25(6):e108-e119.

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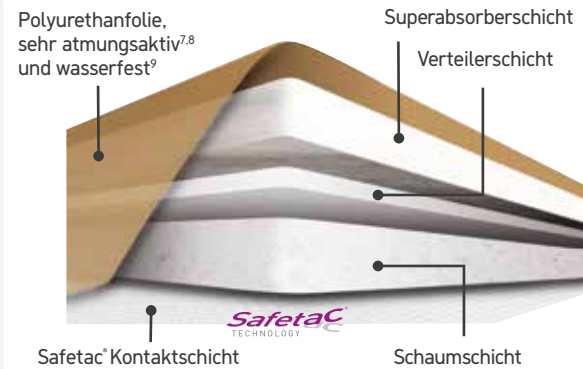
Referenzen:

1. Feili, F. A laboratory evaluation of the fluid management capacities of absorbent foam dressings. Poster presentation at the European Wound Management Association Conference, Vienna, Austria, 2012. 2. Woo, K., Coutts, P.M., Price, P., Harding, K., Sibbald, R.G. A randomised crossover investigation of pain at dressing change comparing two foam dressings *Advances in Skin and Wound Care* 2009;22(7):304-310. 3. Meaume S. et al. A study to compare a new self-adherent soft silicone dressing with a self-adherent polymer dressing in stage II pressure ulcers. *Ostomy Wound Management*, 2003. 4. Curtis K, Bennett P, Heywood E, Marshall L. Peripherally inserted central catheter cushioning: a pilot study comparing gauze with silicone foam. *Clinical Journal of Oncology Nursing* 2015; 19(3):253-256, 2015. 5. Santamaria, N., Gerditz, M., Kapp, S., Wilson, L., Gefen, A. A randomised controlled trial of the clinical effectiveness of multi-layer silicone foam dressings for the prevention of pressure injuries in high-risk aged care residents: The Border III Trial. *International Wound Journal* 2018 DOI: 10.1111/iwj.12891. 6. Kalowes, P., Messina, V., Li, M. Five-layered soft silicone foam dressing to prevent pressure ulcers in the intensive care unit. *American Journal of Critical Care* 2016;25(6):e108-e119. 7. Mölnlycke Health Care data on file report LIMS 20190121-007 (T-1091, Absorption Capacity; T-1099, Joint strength; T-1068, Fluid handling capacity). 8. Mölnlycke Health Care laboratory report 2016-0241. Data on file. 2016. 2018. 9. T-1083 Waterproofness, Mölnlycke Health Care lab report. 20180613-001.

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Referenzen:

1. Morris C, Emsley P, Marland E, Meuleneire F, White R. Use of wound dressings with soft silicone adhesive technology. Paediatric Nursing, April 2009, Vol21, No3:38-43 2. White R., A Multinational survey of the assessment of pain when removing dressings. Wounds UK 2008; Vol 4, No 1, 1-6 3. Randomised Case Study to Evaluate a New Soft Silicone Dressing, Mepilex® Border Lite, for Patients with Diabetic Foot Ulcer. Clinical investigation report 4. Product Manual - Conformability PD-570116 T-1086 Conformability, EN 13726- 5.External Report: FE simulation: Comparison between MxB Flex Lite and Allevyn Gentle Lite, PD-567621 Data on file 6. Product Manual - Impact of Spreading Layer, PD-570117. Data on file. Data on file 7. Product Manual - Retention, PD-587750. Data on file. 8. Product Manual - Fluid Handling Capacity, PD-570308. Data on file. 9.Product Manual - MxB Flex Lite_EM Showerproof PD-580483r ev 01

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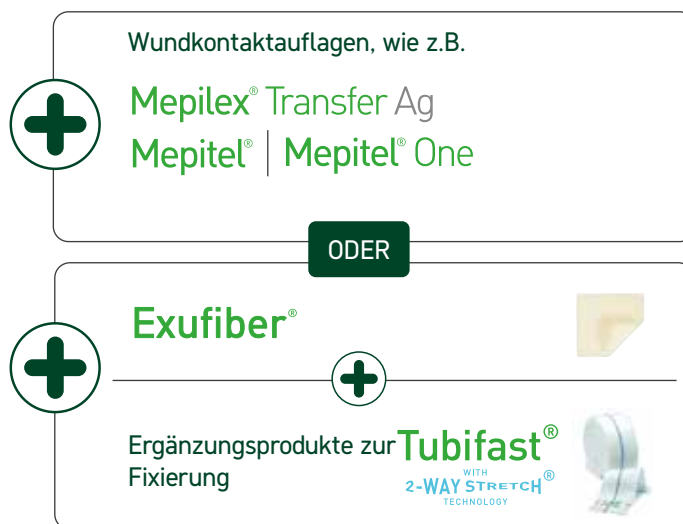
Referenzen:
1. Retrospective Registry Study Report Mesalt, Melgisorb Plus and Mextra, NetHealth. 2019, PD-561171 2. Chadwick P. Haycocks S. Mextra Superabsorbent: case Study compendium - A clinical evaluation of new Superabsorbent dressing in the treatment of a range of foot and leg ulcers. Case reports - Mölnlycke Health Care. 2013. 3. Mextra Superabsorbent Survey 2016, PD-522187 4. Mextra Superabsorbent Survey, 2014, PMS Report PD 486954 (Raw data PD-486953). 5. Product Manual, Mextra Superabsorbent - Absorption and retention of exudate containing bacteria. PD-592501 6. Bacteria trapping testings R-20210212-03 7. Wiegand, C., Hipler, U. A superabsorbent polymer-containing wound dressing efficiently sequesters MMPs and inhibits collagenase activity in vitro. J Mater Sci: Mater Med 24, 2473-2478 (2013) doi:10.1007/s10856-013-4990

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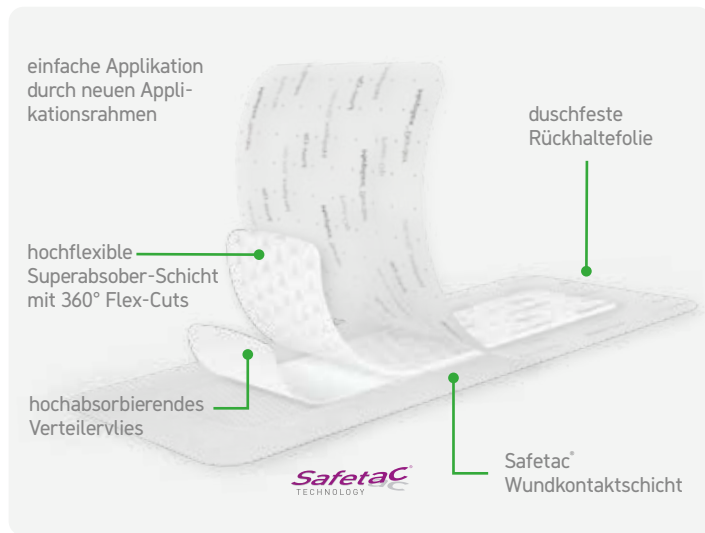
Referenzen:

1. Incision care and dressing selection in surgical wounds: Findings from a series of international meetings. Sandy-Hodgetts K., Morgan-Jones R. Wounds International, 2022. 2. Van Overschelde P, et al. Undisturbed wound healing: a single centre retrospective study investigating patient reported outcomes and clinical validity of extended dressing wear time for incisional healing following orthopaedic surgery, the ARCTIS study. Journal of Wound Care, 2024. 3. Beele H. et al. A prospective randomized controlled clinical investigation comparing two post-operative wound dressings used after elective hip and knee replacement; Mepilex® Border Post-Op versus Aquacel® Surgical. International Journal of Orthopaedic and Trauma Nursing, 2020. 4. Zarghooni K. et al. Is the use of modern versus conventional wound dressings warranted after primary knee and hip arthroplasty? Acta Orthopaedica Belgica, 2015. 5. Bredow J. et al. Evaluation of Absorbent Versus Conventional Wound Dressing. A Randomized Controlled Study in Orthopedic Surgery. Deutsche Arzteblatt International, 2018. 6. Southern K et al. Compliance with the ISWCAP checklist in a THR ERAS pathway – A snapshot audit. Clinical

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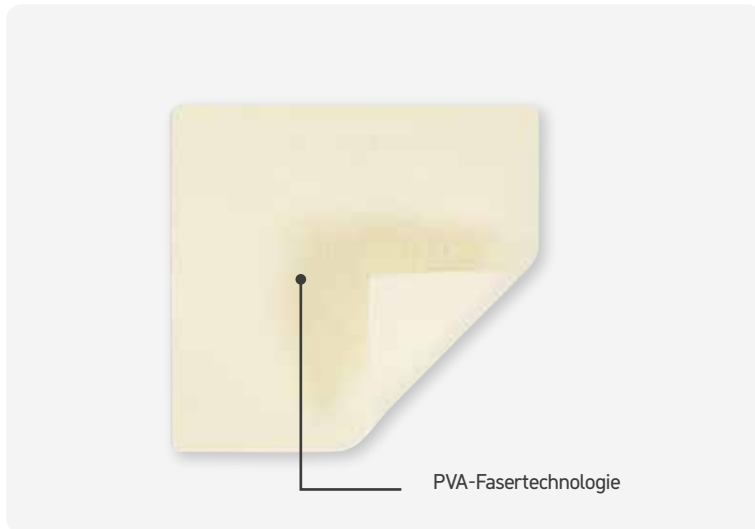
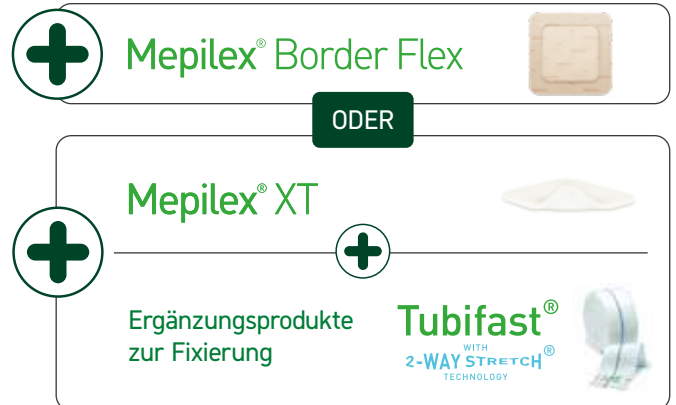
Referenzen:

1. Chadwick P, McCardle J. Open, non-comparative, multicenter post clinical study of the performance and safety of a gelling fibre wound dressing on diabetic foot ulcers. Journal of Wound Care 2016; 25(4): 290-300 2.Smet, S., Beele, H., Saine, L., Suys, E., Henrickx, B. Open, non-comparative, multi centre post market clinical follow-up investigation to evaluate performance and safety on pressure ulcers when using a gelling fibre dressing as intended. Poster Presentation at European Pressure Ulcer Advisory Panel Conference, 2015, Ghent, Belgium 3. Surgical Materials Testing Laboratory. BS EN 13726-1:2002: Test methods for primary wound dressings. Mölnlycke Health Care. Data on file. (2014). 4.Joergensen B, Blaise S, Svensson A-S. A randomised, open-label, parallel-group, multicentre, comparative study to compare the efficacy and safety of Exufiber® with Aquacel® Extra™ dressings in exuding venous and mixed aetiology leg ulcers. Int Wound J. 2022; 19(S1): 22-38 5. Lev-Tov H, Picchietti A. Exufiber Ag+01: A clinical investigation to study the effect of Exufiber Ag+ and other gelling fibre dressings on wound exudate and bioburden in medium to high exuding wounds; Mölnlycke Clinical Investigation Report. PD573989 Rev 01. 2019. Data on file 6. Lustig A, Alves P et al. The sorptivity and durability of gelling fibre dressings tested in a simulated sacral pressure ulcer system. International Wound Journal 2020; doi:10.1111/iwj.13515 [Epub ahead of print] 7. Davies, P., McCarty, S., An in-use product evaluation of a gelling fibre dressing in wound management. E-poster presentation at Wounds UK Conference, 2017, Harrogate, United Kingdom

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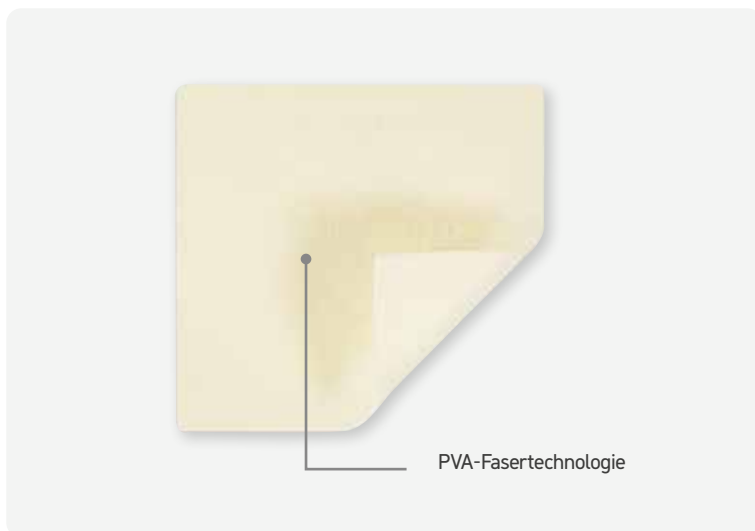
Referenzen:

1. Mölnlycke Health Care. NAMSA reference 10C_39972.01. (Modified) ISO 20743:2007 Determination of antibacterial activity of antibacterial finished products. Data on file, 2010.
2. Mölnlycke Health Care. NAMSA reference 09C_29253.01, 09C_29252.01, 09C_29253.02. ASTM E2149-01 Determining the antimicrobial activity of immobilized antimicrobial agents under dynamic contact conditions. Data on file, 2009.
3. Zurcher S, Krapfl L, Burds A. Reducing postoperative surgical site infections in coronary artery bypass graft patients. Poster presentation at the 45th Annual Conference of the Wound Ostomy and Continence Nurses Society, Seattle, WA, US 2013.
4. Mölnlycke Health Care, Data on File

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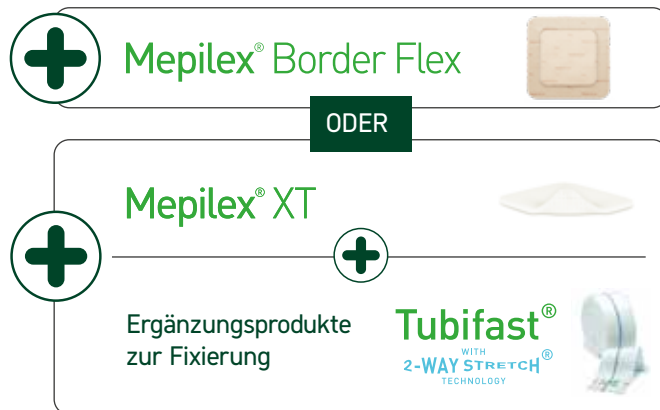
Referenzen:

1. Mölnlycke Health Care, Lev-Tov H. A clinical investigation to study the effect of Exufiber Ag+ and other gelling fibre dressings on wound exudate and bioburden in medium to high exuding wounds. Mölnlycke Data on file, 2018 2. Mölnlycke Health Care, Exufiber Ag+: Physical properties over time (Mölnlycke data on file, 2019) 3. Mölnlycke Health Care, Absorption of thick exudate Exufiber Ag+ vs. Aquacel Ag Extra. (Data on file, 2015.) 4. Mölnlycke Health Care, Exufiber Ag+ Wound Healing Study - Report (Data on file, 2015.) 5. Mölnlycke Health Care, Exufiber Ag+: Competitor tests - Physical (Data on file, 2019) 6. Mölnlycke Health Care, Exufiber Ag+: Challenge Ag: Retention test on Aquacel competitors. (Data on file, 2016) 7. Gil et al. 2017, Evaluation of a Gelling Fiber Dressing with Silver to Eliminate MRSA Biofilm Infections and Enhance the Healing. Poster presented at the Symposium on Advanced Wound Care Spring meeting/Wound Healing Society (WHS) Annual Meeting 2017, Apr 05 - 09, 2017, San Diego, CA, USA. 8. Davis, S. C., Li, J., Gil, J., Head, C., Valdes, J., Glinos, G. D., Solis, M., Higa, A. and Pastar, I. (2019), Preclinical evaluation of a novel silver gelling fiber dressing on Pseudomonas aeruginosa in a porcine wound infection model. Wound Rep Reg. 27: 360-365. 9. Hamberg K., Gerner E. and Falkbring S., Mölnlycke Health Care, Gothenburg, Sweden. "Antimicrobial effect of a new silver-containing gelling fibre dressing against common wound pathogens". Poster presented at the Symposium on Advanced Wound Care Spring meeting/Wound Healing Society (WHS) Annual Meeting 2017, Apr 05 - 09, 2017, San Diego, CA, USA. 10. Hamberg K., Gerner E. and Falkbring S., Mölnlycke Health Care, Gothenburg, Sweden. "In vitro evaluation of the antimicrobial effect of silver-containing fibre dressings". Poster presented at the Symposium on Advanced Wound Care Spring meeting/Wound Healing Society (WHS) Annual Meeting 2017, Apr 05 - 09, 2017, San Diego, CA, USA. 11. Mölnlycke Health Care, CE: Performance of Exufiber Ag+ in vitro: Antimicrobial effect, silver release kinetics and minimal effective concentration. (Data on file, 2016.) 12. Bjarnsholt T, Eberlein T, Malone M, Schultz G. Management of wound biofilm Made Easy. London: Wounds International 2017

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Referenzen:
 1 Schweiger H, Smith D, Cruse CW et al. An open, non-controlled, single-centre, clinical investigation to evaluate efficacy when using a soft silicone wound contact layer containing silver. Poster presentation at the 15th European Burns Association Congress, Vienna, Austria, 2013. 2 Dhatriya K, Gooday C, Franke B, Pilling T, Flanagan A, Zeidan L. (2016): An open, non-comparative, multicentre evaluation of performance and safety using an antimicrobial exudate transfer dressing on diabetic foot ulcers: a case series. In Journal of wound care 25(5), pp. 256-265. 3 Mölnlycke Health Care, Data on file, 2016. Arenberger P. A prospective open noncontrolled clinical investigation to evaluate the performance and safety of an antimicrobial foam transferring dressing in patients with malignant wounds. Clinical Investigation Report MxT Ag 03. 2016. 4 Scheske J, Würfel T, Kudrina M. GVV Registry Study Mepilex Transfer Ag. 2016. 5 T-2133 Passage of fluid through wound contact layers. Report no.: 20190306-003 6 T-1112 Antimicrobial effect of Mepilex Transfer Ag against P. aeruginosa, S. aureus and C. albicans at 30 minutes, 3 hours, 24 hours and 8 days of incubation time. Report no.: 20111111-003 7 T-1112 Antimicrobial effect of Mepilex Transfer Ag against 16 strains at 30 minutes, 3 hours, 24 hours and 8 days of incubation time, performed by NAMSA. Report no.: 20120206-005 8 T-1112 Re-test of antifungal effect of Mepilex Transfer Ag against 2 yeast strains at 30 minutes, 3 hours, 24 hours and 8 days of incubation time, performed by NAMSA. Report no.: 20120323-001 9 T-1116 Antimicrobial effect of Mepilex Transfer Ag in combination with secondary dressing against P. aeruginosa, S. aureus and C. albicans at 24 hours. Report no.: 20120116-005 10 T-1116 Antimicrobial effect of Mepilex Transfer Ag in combination with secondary dressing against P. aeruginosa, S. aureus and C. albicans at 24 hours. Report no.: 20120207-010 11 T-1112 Antimicrobial effect of Mepilex Transfer Ag against P. aeruginosa, S. aureus and C. albicans at 24 hours of incubation time. Report no.: 20110328-006

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Mextra®
Superabsorbent

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Tubifast®
WITH 2-WAY STRETCH TECHNOLOGY

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Referenzen:
1. Chadwick P, Taherinejad F, Hamberg K, Waring M. Clinical and scientific data on a silver-containing soft-silicone foam dressing: an overview. Journal of Wound Care 2009;18(11):483-491. 2. Mölnlycke Health Care. Mepilex Ag: Determining antimicrobial activity of immobilized antimicrobial agents (ASTM-E2149-01) (NAMS). Data on file, 2006. 3. Mölnlycke Health Care. Determination of sustained antimicrobial effect of Mepilex Ag using CZOI (Clinical Bacteriology). Data on file, 2006. 4. Mölnlycke Health Care. Sustained antimicrobial effect up to 8 days in ASTM-E2149-01. Data on file, 2007. 5. Mölnlycke Health Care. Antimicrobial activity against a broad range of bacteria and fungi. Data on file, 2007. 6. Tang H, Lv G, Fu J, Niu X, Li Y, Zhang M, et al. An open, parallel, randomized, comparative, multicenter investigation evaluating the efficacy and tolerability of Mepilex Ag versus silver sulfadiazine in the treatment of deep partial-thickness burn injuries. Journal of Trauma and Acute Care Surgery 2015;78(5):1000-1007. 7. Rigo C, Roman M, Munivra I, Vindigni V, Azzena B, Barbante C, et al. Characterization and evaluation of silver release from four different dressings used in burns care. Burns 2012;38(8):1131-1142. 8. Gee Kee EL, Kimble RM, Cuttle L, Khan A, Stockton KA. Randomized controlled trial of three burns dressings for partial thickness burns in children. Burns 2015;41(5):946-955. 9. Silverstein P, Heimbach D, Meites H, Latenser B, Mozingo D, Mullins F, et al. An open, parallel, randomized, comparative, multicenter study to evaluate the cost-effectiveness, performance, tolerance, and safety of a silver-containing soft silicone foam dressing (intervention) vs silver sulfadiazine cream. Journal of Burn Care & Research 2011;32(6):617-626. 10. Tang H, Lv G, Fu J, Niu X, Li Y, Zhang M, et al. An open, parallel, randomized, comparative, multicenter investigation evaluating the efficacy and tolerability of Mepilex Ag versus silver sulfadiazine in the treatment of deep partial-thickness burn injuries. Journal of Trauma and Acute Care Surgery 2015;78(5):1000-1007. 11. Richards K, Chadwick P. Addressing local wound infection with a silver-containing, soft silicone foam dressing: a case series. The Diabetic Foot Journal 2011;14(2):90-95.

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Referenzen:

1. Collin O. Use of Mepitel One dressing following hand surgery: a case study series. Poster presentation at Wounds UK Conference, United Kingdom, 2009. 2. Dahlstrom K.K. A new silicone rubber dressing used as a temporary dressing before delayed split skin grafting: A prospective randomised study. Scandinavian Journal of Plastic and Reconstructive Surgery and Hand Surgery, 1995. 3. Brölmann F., et al. Randomized clinical trial of donor-site wound dressings after split-skin grafting. British Journal of Surgery, 2013. 4. Terrill P.J. A comparison of three primary non-adherent dressings applied to hand surgery wounds. J Wound Care, 2000. 5. O'Donovan DA et al. The role of Mepitel silicone net dressings in the management of fingertip injuries in children. Journal of Hand Surg, 1999. 6. Platt A.J et al. A comparative study of silicone net dressing and paraffin gauze dressing in skin-grafted sites. Burns, 1996. 7. Eid R. et al. SIMP 01 - A prospective, open, non-controlled proof-of-concept investigation to evaluate the function of a new silicone cover film and an on-top suction device to be used in a NPWT system in the treatment of acute and chronic wounds. Mölnlycke Clinical Investigation Report, 2011. 8. Adamietz, I. A. et al. Effect of Self-Adhesive, Silicone-Coated Polyamide Net Dressing on Irradiated Human Skin. Radiation Oncology Investigations, 1995. 9. Bugmann P.H. et al. A silicone coated nylon dressing reduces healing time in burned paediatric patients in comparison with standard sulfadiazine treatment: a prospective randomized trial. Burns, 1998. 10. Campanella SD, et al. A randomised controlled pilot study comparing Mepitel and SurfaSoft on paediatric donor sites treated with Recell. Burns, 2011;37(8):1334-42. 11. Mölnlycke Health Care. Design Verification Report, Mepitel One, PD-557646. Data on file. 12. Patton P. et al. An open, prospective, randomized pilot investigation evaluating pain with the use of a soft silicone wound contact layer vs. a lipidocolloid wound contact layer (UrgoTul) in the treatment of acute wounds. International Wound Journal, 2017. 14. Collin O. Use of Mepitel One dressing following hand surgery: a case study series. Poster presentation at Wounds UK Conference, United Kingdom, 2009.

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Referenzen:
1. Meuleneire F. A vapour-permeable film dressing used on superficial wounds. British Journal of Nursing. 2014. 2. Möller PK, et al. Breast cancer patients report reduced sensitivity and pain using a barrier film during radiotherapy – A Danish intra-patient randomized multicentre study. Technical Innovations and Patient Support in Radiation Oncology. 2018. 3. Herst P, et al. Prophylactic use of Mepitel Film prevents radiation-induced moist desquamation in an intra-patient randomised controlled clinical trial of 78 breast cancer patients. Radiotherapy and Oncology. 2014 4. Mölnlycke Health Care. Data on file. 5. Wooding H, et al. The effect of Mepitel Film on acute radiation-induced skin reactions in head and neck cancer patients: A feasibility study. British Journal of Radiology. 2018.

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Referenzen:
1. Van Overschelde P, et al. Undisturbed wound healing: a single centre retrospective study investigating patient reported outcomes and clinical validity of extended dressing wear time for incisional healing following orthopaedic surgery, the ARCTIS study. Journal of Wound Care, 2024 2 Internal laboratory test (T-1191) reported in Design Validation Record. Data on file. 3 Product System Description. Data on file. 4 Internal laboratory test (T-1180) reported in Design Verification Report Data on file. 5 Van Overschelde, P. et al. A randomised controlled trial comparing two wound dressings used after elective hip and knee arthroplasty. Poster presentation at 5th Congress of the WUWHS, Florence, Italy, 2016. 6 Silverstein P. et al. An open, parallel, randomized, comparative, multicenter study to evaluate the costeffectiveness, performance, tolerance and safety of a silver-containing soft silicone foam. Journal of Burn Care and Research, 2011. 7 Gee Kee E.L. et al. Randomized controlled trial of three burns dressings for partial thickness burns in children. Burns, 2014. 8 David F. et al. A randomised, controlled, non-inferiority trial comparing the performance of a soft silicone-coated wound contact layer (Mepitel One) with a lipidocolloid wound contact layer (UrgoTul) in the treatment of acute wounds. International Wound Journal, 2017. 9 Elhage KG, Awad ME, Irfan FB, Lumbley J, Mostafa G, Saleh KJ. Closed-incision negative pressure therapy at 125 mmHg significantly reduces surgical site complications following total hip and knee arthroplasties: A stratified meta-analysis of randomized controlled trials. Health Sci Rep. 2022;5:e425. doi:10.1002/hsr.2.425

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Referenzen:

1. Edwards-Jones V, Flanagan M, and Wolcott R. Technological advancements in the fight against antimicrobial resistance. Wounds Int 2015;6(2):47-51. 2. Granudacyn®- Journal of Wound Care Volume 29 No 10 (Suppl 2) October 2020. 3. Severing AL, Rembe JD, Koester V, Stuermer EK. Safety and efficacy profiles of different commercial sodium hypochlorite/hypochlorous acid solutions (NaClO/HClO): antimicrobial efficacy, cytotoxic impact and physicochemical parameters in vitro Journal of Antimicrobial Chemotherapy, Volume 74, Issue 2, February 2019, Pages 365-372. 4. RDS 005-0414, ZLG, Anleitung für die Festlegung von Mindestkriterien zur mikrobiologischen Reinheit von Medizinprodukten. 5. Efficacy and safety of neutral pH superoxidised solution in severe diabetic foot infections - Fermín R Martínez-De Jesu's, Antonio Ramos-De la Medina, José María Remes-Troche, David G Armstrong, Stephanie C Wu, Jose Luis Lazaro Martinez, Juan V Beneit-Montesinos - International Wound Journal.

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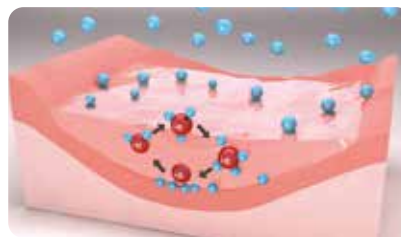
Referenzen:
1. Ingo, S., Petri, M., Dissemond, J., Joachim, K. Photoacoustic tomography for non-invasive diagnostic in patients with chronic venous leg ulcers treated with hemoglobin spray – first results of a prospective clinical study. Oral presentation, EWMA 2015, London. 2. Dissemond, J., Kröger, K., Storck, M., Risse, A., Engels, P. Topical oxygen wound therapies for chronic wounds: a review. Journal of Wound Care, 2015, Feb, 24(2);53-63. 3. Gottrup F, Dissemond J, Baines et al. Use of oxygen therapies in wound healing, with special focus on topical and hyperbaric oxygen treatment. J Wound Care, 2017; 26(5), Suppl, S1-S42. 4. Hunt, SD, Elg, F. Clinical effectiveness of hemoglobin spray (Granulox) as adjunctive therapy in the treatment of chronic diabetic foot ulcers. Diabetic Foot & Ankle, 2016; 7(1):33101. https://doi.org/10.3402/dfa.v7.33101. 5. Hunt, SD, Elf, F, Percival, SL. Assessment of clinical effectiveness of haemoglobin spray as adjunctive therapy in the treatment of sloughy wounds. J Wound Care, 2018; 27(4): 210-219. 6. Elg, F, Bothma, G. Cost effectiveness of adjunct haemoglobin spray in the treatment of hard-to-heal wounds in a UK NHS primary care setting. J Wound Care, 2019; 28(12):5. 7. Hunt, SD, Elg, F. The clinical effectiveness of haemoglobin spray as adjunctive therapy in the treatment of chronic wounds. J Wound Care, 2017; 26(9). 8. Brüggengjürgen, B, Hunt, SD, Eberlein, T. Wound management in diabetic foot ulcer (DFU) – incremental cost-analysis of treating diabetic neuropathic foot lesions with adjunct hemoglobin contact spray in Germany. Gesundheitsökonomie und Qualitätsmanagement, 2017; 22: 1-8. 9. Petri, M., Stoffels, I., Griewank, K., Jose, J., Engels, P., Schulz, A., Pötzschke, H., Jansen, P., Schadendorf, D., Dissemond, J., Klode, J. Oxygenation Status in Chronic Leg Ulcer After Topical Hemoglobin Application May Act as a Surrogate Marker to Find the Best Treatment Strategy and to Avoid Ineffective Conservative Long-term Therapy. World Molecular Imaging Society, 2017. 10. Scholander, PF. Oxygen transport through hemoglobin solutions. Science, 1960 Feb 26;131(3400):585-90. 11. R. Strohal et al. Expertenkonsens zu praxisrelevanten Aspekten der Wundtherapie mit einem Hämoglobin-Spray. Wundmanagement 2016; 5: 276-284.

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Referenz zu Mepore Pro: 1. Mölnlycke Health Care; Data on file, 2019

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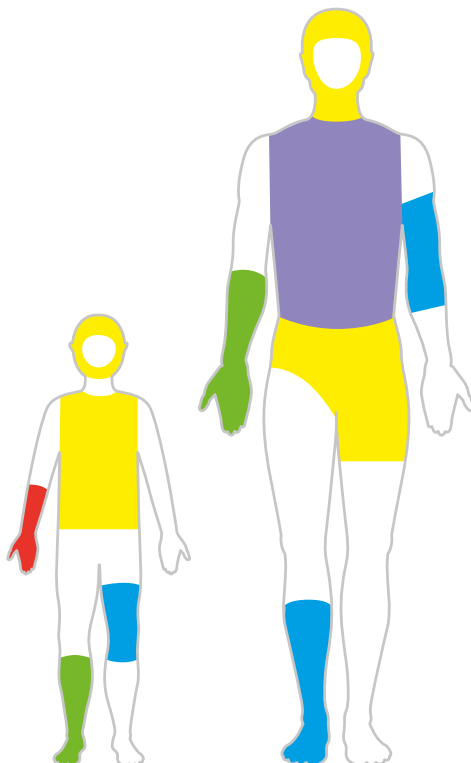
Referenzen zu Tubifast:

1. Eytier C, Gazeau E, Beneteau G, Verfaillie G. Convenience and tolerance of the combination of a soft silicone foam dressing and a two-way stretch tubular bandage in the management of local wounds. Journal des plaies et cicatrisations 2013;18(88):38-44. 2. Mölnlycke Health Care. Data on file.

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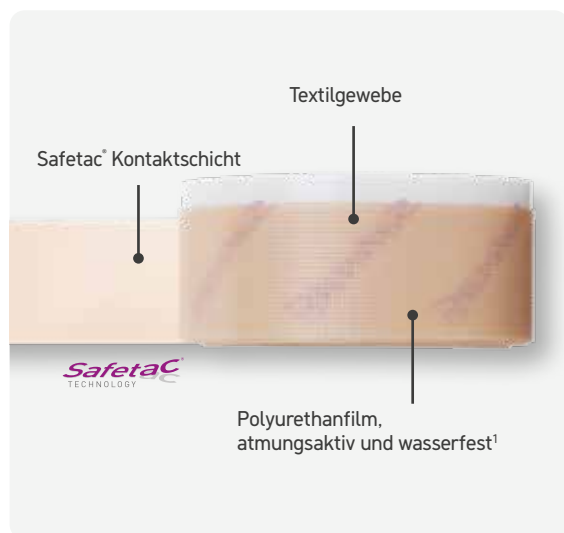
Referenzen:
1. Mölnlycke data on file. 2. Livesley J, Richardson S. Securing methods for peripheral cannulae. Nursing standard (Royal College of Nursing (Great Britain)) : 1987. 1993;7(31):31-4. 3. Kanapathy M, Mosahebi A. Comparative study on the donor site aesthetic outcome between epidermal graft and split-thickness skin graft. International Wound Journal. 2019;16(2):354-9

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Referenz:

1. White, R. A multinational survey of the assessment of pain when removing dressings. Wounds UK 2008; 4: 1, 14-22 2. Colom Majan, J.L. Evaluation of a self-adherent soft silicone dressing for the treatment of hypertrophic post-operative scars. Journal of Wound Care 2006;15(5):193-196 3. Mustoe, T.A., Cooter, R.D., Gold, M.H., Hobbs, R., Ramelet, A.A., Shakespeare, P.G. et al., Stella, M., Teot, L., Wood, F.M., Ziegler, U.E., for the International Advisory Panel on Scar Management. International clinical recommendations on scar management. Plastic and Reconstructive Surgery 2002;110 (2): 560-571. 4. Colom Majan, J.L. Evaluation of a self-adherent soft silicone dressing for the treatment of hypertrophic post-operative scars. Journal of Wound Care 2006;15(5):193-196. 5. Meuleneire, F. Using a soft silicone dressing (Mepiform) to prevent scarring in an acute traumatic wound. Wounds UK 2007;3(4):134- 135. 6. Saulsberry, C.M., Abney, S.B., Landry, K.N. An evaluation of a new silicone dressing for the management of hypertrophic scars. Poster presentation at the 12th Symposium on Advanced Wound Care and Medical Research Forum on Wound Repair, Anaheim, California, United States of America, 1999. 7. Data on file, Mölnlycke Health Care, 2015.

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Referenz:

1. Mölnlycke Health Care. Data on file. 2019. 2. Apelqvist J, Ragnarson Tennvall G.; Larsson J. Topical treatment of diabetic foot ulcers. An economic analysis of treatment alternatives and strategies. Diabetic medicine. 1995, 12(2):123-128. 3. Apelqvist J, Ragnarson-Tennvall G. Cavity foot ulcers in diabetic patients. A comparative study of cadexomer iodine ointment and standard treatment. An economic analysis alongside a clinical trial. Acta Derm Venereol. 1996, 76(3): 231-235. 4. Brown-Etris M, Myers RB, Pasceri P. A new generation of gauze dressings. Ostomy/Wound Management, The Journal for Extended Patient Care Management. 1991;34. 5. Dinstl K, Platthy G, Zita L. Case study: Suitability of Hypergel, Mesalt and Normigel in the treatment of leg ulcers. Study report 1995. 6. Hirsch DJ, Jindal KK. Local care of Staphylococcus aureus exit-site infection precludes antibiotic use. Peritoneal Dialysis International. 2003;23(3):301-2. 7. Huang MR, Huang LE, Li MY, Guo SY. Nursing of Postoperative Infected Wound of Children with Sacrococcygeal Teratoma. Journal of Nursing. 2011, 16:14. 8. Khanduja KS, Padmanabhan A, Kerner BA, Wise WE, Aguilar PS. Reconstruction of rectovaginal fistula with sphincter disruption by combining rectal mucosal advancement flap and anal sphincteroplasty. Diseases of the colon & rectum. 1999, 42(11):1432-1437. 9. Liu Z, Xiong Z, Wu J, Wang F. Clinical curative effect of Mesalt combined with Mepilex dressing in postoperative infection of inguinal hernia. Med Sci Monit. 2015 Apr 9;21:1038-42. 10. Qin YU. Dressing Effect of Mesalt Dressing on Infected Wounds. Journal of Kunming Medical University. 2012, 10:54. 11. Qiuping Z, Wei Z. Application of mesalt dressing in treating infected wound. Journal of Nursing Science. 2010, 22:36. 12. Valente MA, Khanduja KS. Layered surgical repair of traumatic cloacal deformities. Technical details and functional outcomes. Techniques in coloproctology. 2012, 16(2):153-156. 13. Zhou S, Ling W, Xia D. Observation of the effect of Urgotul parcel technology used for treatment of surface abscess incision and drainage of the wound. Chinese Journal of Practical Nursing. 2018;34(3):187-91. 14. Mölnlycke Health Care. Data on file. 2020

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Referenz:

1. Mölnlycke Health Care. Data on file. 2019 2. Mölnlycke Health Care. Data on file. 1998

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Die Namen Mölnlycke Health Care und das Safetac Logo sowie alle genannten Produkte sind weltweit eingetragene Marken eines oder mehrerer Mitglieder der Mölnlycke Health Care Unternehmensgruppe. ©2023 Mölnlycke Health Care. Alle Rechte vorbehalten. DWC0128 / Stand: September 2024